

Klarity[®] BiteLok[™]

Higher precision for treatment accuracy




RNB-100

Klarity BiteLok[™]

The Klarity BiteLok[™] is an improved, innovative design for increased head and neck immobilization. The BiteLok[™] is made of a rigid, breathable mouthpiece. It indexes to the patient's mouth using fast-setting dental putty for a secure and comfortable fit. No water or warming is needed to activate the putty. Once the BiteLok[™] is indexed to the patient's mouth, a thermoplastic mask is placed overtop. While the mask is still warm, a hands-free Quickfit tool forms the mask around the edges of the BiteLok[™] mouthpiece. The result is a secure, reproducible mouthpiece that indexes to both patient and mask, that can also be removed for cleaning.

The Klarity BiteLok[™] features a smaller profile than any other mouthpiece, and the wide airway passage allows for better breathability. The BiteLok[™] is a comfortable size for most patients, even those with smaller mouths. The BiteLok[™] is a significant advancement in comfortable, accurate positioning for head and neck treatments.

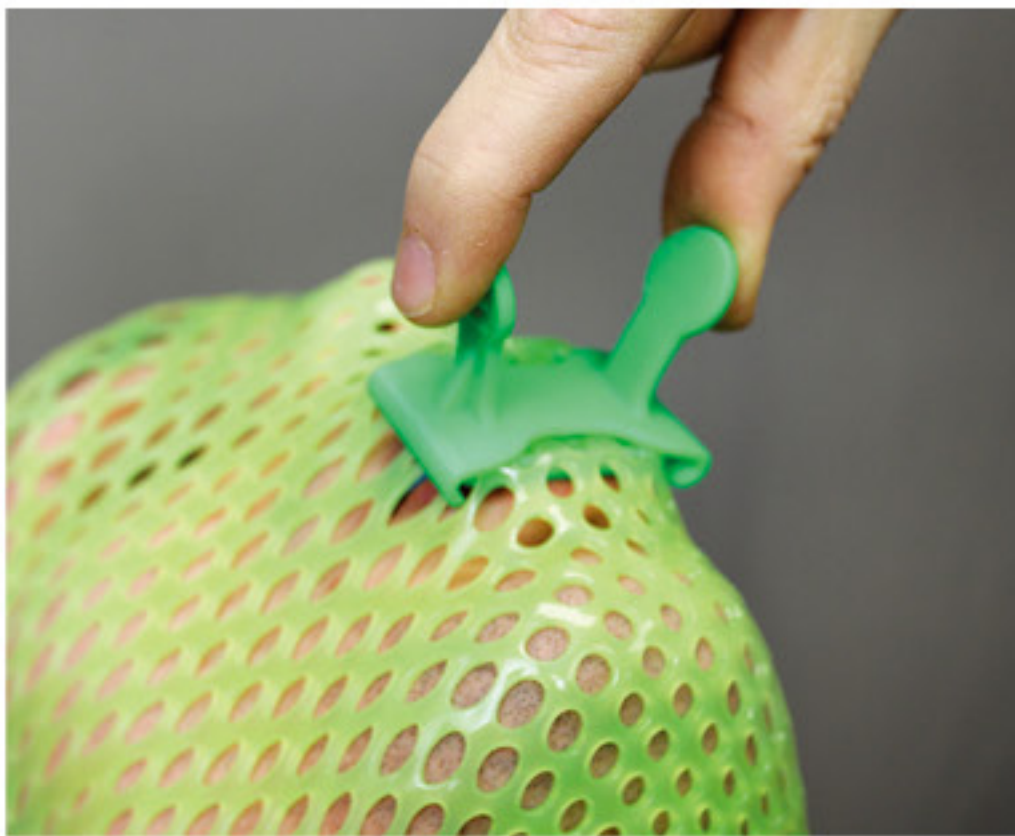

QuickFit tool indexing a thermoplastic mask over the BiteLok[™]



OPTIONAL INDEXING CHOICES FOR CUSTOMIZED POSITIONING

The BiteLok™ can be equipped with a tongue depressor or tongue diverter for additional mouth placement options.

The BiteLok™ Tongue Depressors come in two different sizes, both offering a comfortable fit that keeps the tongue in a down position without obstructing the airway. The tongue depressors securely snaps into the BiteLok™. A removal tool is included if the tongue depressors needs to be removed or changed.


The BiteLok™ Tongue Diverter is an additional option for a customized fit. Again, the tongue diverter snaps easily into the BiteLok™, and can also be removed with the removal tool.

- Easily breathable mouthpiece
- Indexes to mask
- Removable for easy cleaning
- Smaller profile for better patient comfort
- Optional Tongue Depressor and Tongue Diverter

OPTIONAL COMPONENTS:

Each Klarity BiteLok™ comes as a complete set with the following optional components:


Tongue Depressor

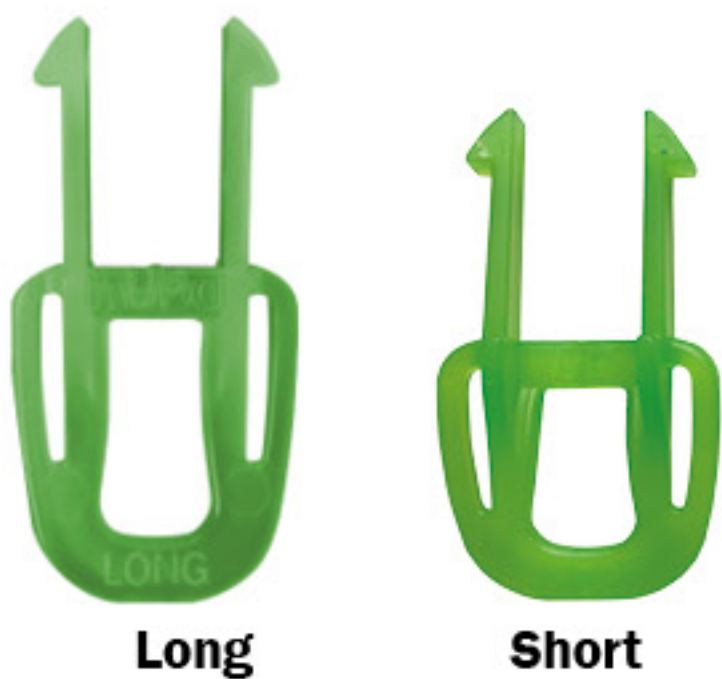

Tongue Depressor



Tongue Diverter

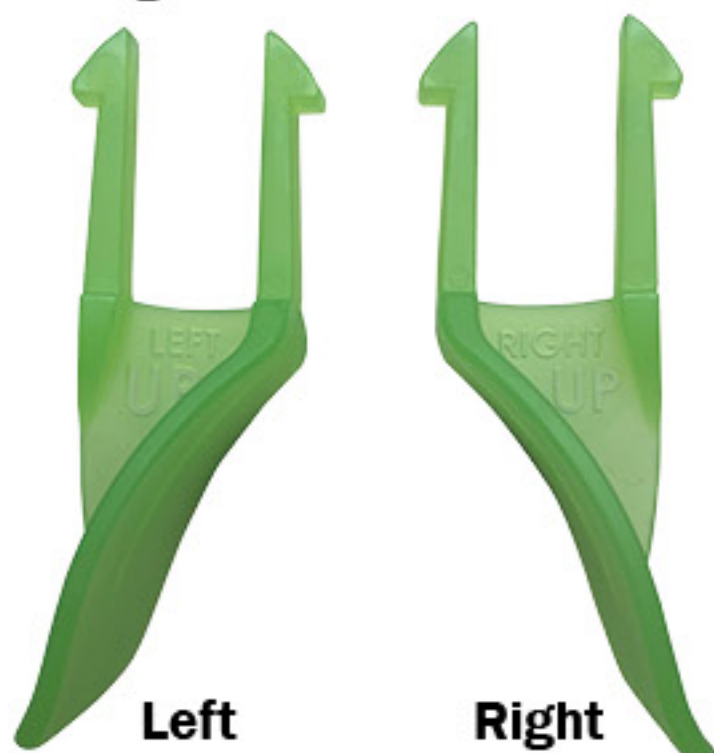

Tongue Diverter

