

Klarity® **ATHENA**™

Supine Breast Solutions

- Enables better personalized immobilization
- More reliable with max load limit 200 kgs
- MR compatible version available



Klarity ATHENA



Klarity ATHENA™ Supine Breast Solutions

- Ideal for the combination use with masks, vacuum bags or Accucushions for better personalized immobilization.
- Enables patient positioning head sideways to the contralateral side.
- Maximum load limit up to 200kgs.

Klarity ATHENA™ Breast System with Breastboard + Wing Board

Features:

- The Wing Board can be easily and steadily fixed onto the Breastboard, providing comfortable overhead arm positioning.
- The System can use with S-type head mask, head & thorax mask and a dedicated vacuum bag.
- Elevation support with angles at 5°, 10°, 15° and 20°.

* The System doesnot include Arm Supports and Hand-grips.



Leg Support Cushion



Klarity ATHENA™ Breastboard

Features:

- The Breastboard is equipped with easy-to-fix and easy-to-adjust Arm Supports and Hand-grips, providing comfortable one-arm or two-arm overhead positioning.
- A special headrest for patient positioning head sideways to the contralateral side is included.
- Elevation support with angles at 5°, 10°, 15° and 20°.



Klarity ATHENA™ Wing Board




Code: R632-NTCF
Material: Carbon fiber



Code: R632-NTKF
Material: Kevlar




Klarity ATHENA™ Breast System

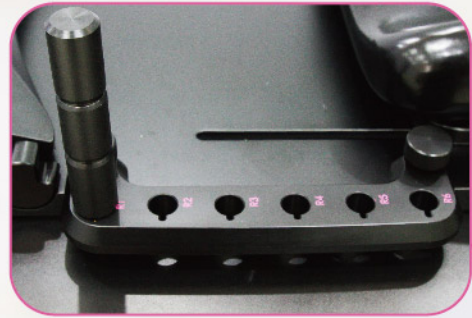

Adjustable Hand-grips at the bottom provide comfortable positioning for patients of different heights and enable one arm overhead positioning.



The T-Handle of the Wing Board is easy to adjust along the Y and Z axis.

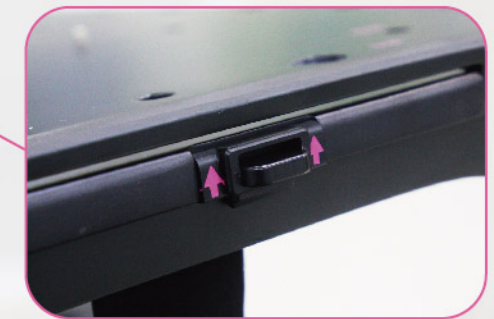

The Elevation Supports are easy to be secured and unlocked.



A special adapter is designed for easy mask/cushion fixation.



Zero-pressure in hip thanks to the V-shaped structure greatly improves patient comfort and positioning stability.



The Klarity ATHENA™ is well accommodated in a 70cm bore when elevated at 5°.





Klarity ATHENA™ Breast System
for Breast & Thorax
Personalized Immobilization
with the use of mask and vacuum bag



Breast & Thorax Personalized Immobilization

Klarity ATHENA™ Breast System



Klarity ATHENA™
+ S-type Head Mask
+ Vacuum Bag

* Ideal for patient positioning head
sideways to the contralateral side

Breast & Thorax Personalized Immobilization

Klarity ATHENA™ Breast System



Klarity ATHENA™
+ Head & Thorax Mask (open at treated region)
+ Vacuum Bag

* Ideal for patient positioning head
sideways to the contralateral side

Compatible Consumables



• Vacuum Bag:
MR7677-25N




• Head & Thorax Mask:
R322-S424C1-L (White, Left side)
R322-S424C1-R (White, Right side)



• S-type Head Mask:
R460SHT (White, 2.4mm)
R461SHT (White, 3.2mm)
RLG460SHT (LiteGreen™, 2.4mm)
RLG461SHT (LiteGreen™, 3.2mm)
RB461-3210SHT (Blue, 3.2mm)

Details

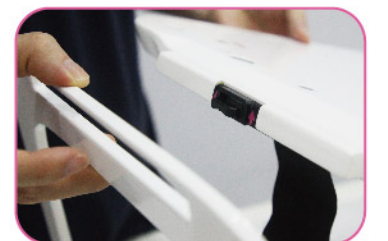

Klarity ATHENA™ Breastboard



The soft pad of the easy-to-fix and easy-to-adjust Arm Support is replaceable and easy to clean.



A special headrest is designed for patient positioning head sideways to the contralateral side.



Zero-pressure in hip
thanks to the V-shaped structure
greatly improves patient comfort and
positioning stability.



Adjustable Hand-grips at the bottom
provide comfortable positioning for
patients of different heights and
enable one arm overhead positioning.



Klarity

Klarity ATHENA™ Breastboard

for Breast & Thorax

Personalized Immobilization

with the use of mask and Accucushion®



Two-arm Overhead Positioning Solution

Klarity ATHENA™ Breastboard

- Klarity ATHENA™
- + Head & Thorax Mask
- + Head & Shoulder Accucushion®
- + Headrest

* Ideal for patient positioning head sideways to the contralateral side



Compatible Consumables




• Head & Shoulder AccuCushion®:
RG560-SH1



• Head & Thorax Mask:
R322-S424C1-L (White, Left side)
R322-S424C1-R (White, Right side)



• S-type Head Mask:
R460SHT (White, 2.4mm)
R461SHT (White, 3.2mm)
RLG460SHT (LiteGreen™, 2.4mm)
RLG461SHT (LiteGreen™, 3.2mm)
RB461-3210SHT (Blue, 3.2mm)

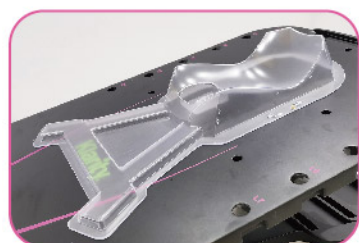

• Dedicated Headrest for AccuCushion®:
R507-T

One-arm Overhead Positioning Solution

Klarity ATHENA™ Breastboard

- Klarity ATHENA™
- + Head & Shoulder Accucushion®
- + Headrest



Details





BETTER PRODUCTS FOR BETTER CARE



Klarity Medical & Equipment (GZ) Co., Ltd.
No.14, 3rd Street Shawan,
GETDD Guangzhou 510730,China

Tel: +86-20-8222 6380

Fax: +86-20-8208 8552

Sales@klarity-medical.com

Version No.: KL-CY-20220720-A0



LinkedIn



YouTube